

Product data sheet






Information

Designation Reference Taupe and chrome grey foldaway shower seat - Large size, height: 550 mm 047625

Product advantages

- + Maximum safety: with automatic stand.
- + Compact in raised position.
- + Easy fixing
- + Corrosion resistant warranty.

Description

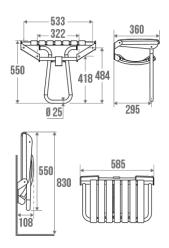

Large size. 7 taupe (Pantone 7591U) polypropylene slats. Grey epoxy-coated aluminium tube. Enables users to sit down safely and comfortably while showering.

Specific height 550 mm, for private use only. Public facilities standard : 500 mm.

Technical characteristics

Finish	Epoxy-coated
Color	Chrome Grey & Taupe
Dimensions	360 x 580 x 550 mm
Diameter Ø	25 mm
Material	Aluminium & Polypropylene
Gross weight	4.700 kg
Packaging	Cardboard box
Maximum Weight Capacity	150 kg

PELLET ASC - ZI DE CHAPOTIN - 240 RUE DES FRÈRES VOISIN - CS 270048 - 69967 CHAPONNAY CEDEX - FRANCE - TEL : +33 (0)4 78 96 82 20 - FAX : +33 (0)4 78 96 70 52

© PELLET ASC – Non contractual document. Please read our General Sales Conditions. Dimensions in mm. All dimensions and specifications are given for information purposes only and can be changed without notice in order to improve our products.



Laboratory tests Made in France Maintenance Guarantee Gencod Custom code Technical sheet Spare parts 225 kg yes yes 10 years 3563390476253 7615200000 S00001530_IND_02 Ref. 047628, Taupe (Pantone 7591U) seat slat. Ref. 047644,Foot cap for shower seat. Ref. 047623, Support slat.

Fixing and maintenance advices

Careful consideration and selection of the appropriate wall fixings according to the wall-mounted support is required to ensure the long life of your accessory. The use of stainless steel screws (Ø 4 mm) is highly recommended to prevent corrosion in damp environment. Before drilling, physically hold the product against the wall at desired position and mark the hole locations. All the fixing points must be respected. After installation and before use, make sure that all fixing screws are tightened securely. The shower seat provides support and added safety if properly fitted in accordance with our technical instructions. For use in seated position. To guarantee the user safety, DO NOT rise up on the seat. When the seat is raised, DO NOT use it as a grab bar. To avoid any risks of pinching out, DO NOT put your fingers in the mobile parts. The product does not contains any dangerous substance for the health. We disclaim all liability for any personal injury resulting from inappropriate installation or use of this seat.

Fiche produit mise à jour le : 10/01/2024