

Product data sheet





Information


Designation White and chrome grey foldaway shower seat

Reference 047631

Product advantages

- + Maximum safety: with automatic stand.
- + Compact in raised position.
- + Easy fixing
- + Corrosion resistant warranty.

Description

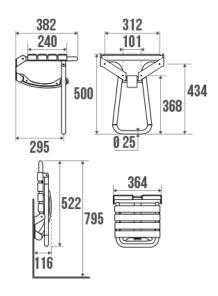

Standard size.

4 white polypropylene slats.

Grey epoxy-coated aluminium tube.

Enables users to sit down safely and comfortably while showering.

Standard height 500 mm.

Technical characteristics

Finish

Color

Dimensions

Diameter Ø Material

Gross weight

Epoxy-coated


Chrome Grey & White

380 x 355 x 500 mm

25 mm

Aluminium & Polypropylene

2.700 kg

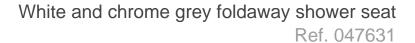



Packaging Four-color process cardboard

 CE Marking : safe load
 150 kg

 Laboratory tests
 225 kg

 Made in France
 yes

 Maintenance
 yes

 Guarantee
 10 years

 Net weight
 2.174 kg

 Gencod
 3563390476314

 Custom code
 7615200000

 Technical sheet
 A13470022 IND 21

Spare parts Ref. 047636: pads for shower seat. Ref. 047638: White seat slat. Ref. 047639: Support slat.

Fixing and maintenance advices

Careful consideration and selection of the appropriate wall fixings according to the wall-mounted support is required to ensure the long life of your accessory. The use of stainless steel screws (Ø 4 mm) is highly recommended to prevent corrosion in damp environment. Before drilling, physically hold the product against the wall at desired position and mark the hole locations. All the fixing points must be respected. After installation and before use, make sure that all fixing screws are tightened securely. The shower seat provides support and added safety if properly fitted in accordance with our technical instructions. For use in seated position. To guarantee the user safety, DO NOT rise up on the seat. When the seat is raised, DO NOT use it as a grab bar. To avoid any risks of pinching out, DO NOT put your fingers in the mobile parts. The product does not contains any dangerous substance for the health. We disclaim all liability for any personal injury resulting from inappropriate installation or use of this seat

Fiche produit mise à jour le : 10/01/2024