

## Product data sheet




### Information

Designation	Grey and chrome grey foldaway shower seat - Large size
Reference	047643

### Product advantages

- + Maximum safety: with automatic stand.
- + Compact in raised position.
- + Easy fixing
- + Corrosion resistant warranty.

### Description

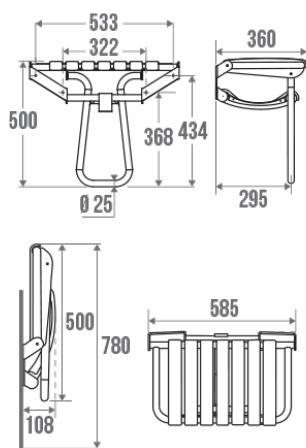

Large size.  
7 Cool Gray (Pantone 7U) polypropylene slats.  
Grey epoxy-coated aluminium tube.  
Enables users to sit down safely and comfortably while showering.

Standard height 500 mm.

### Technical characteristics

Finish	Epoxy-coated
Color	Chrome Grey & Grey
Dimensions	360 x 580 x 500 mm
Diameter Ø	25 mm
Material	Aluminium & Polypropylene
Gross weight	3.600 kg
Packaging	Cardboard box + Label
CE Marking : safe load	150 kg

Laboratory tests	225 kg
Made in France	yes
Maintenance	yes
Guarantee	10 years
Net weight	3.000 kg
Gencod	3563390476437
Custom code	7615200000
Technical sheet	<a href="#">A13600811_IND_18</a>
Spare parts	Ref. 047648, Grey (Pantone Cool Gray 7U) seat slat. Ref. 047644, Foot cap for shower seat.

## Fixing and maintenance advices

---

Careful consideration and selection of the appropriate wall fixings according to the wall-mounted support is required to ensure the long life of your accessory. The use of stainless steel screws (Ø 4 mm) is highly recommended to prevent corrosion in damp environment. Before drilling, physically hold the product against the wall at desired position and mark the hole locations. All the fixing points must be respected. After installation and before use, make sure that all fixing screws are tightened securely. The shower seat provides support and added safety if properly fitted in accordance with our technical instructions. For use in seated position. To guarantee the user safety, DO NOT rise up on the seat. When the seat is raised, DO NOT use it as a grab bar. To avoid any risks of pinching out, DO NOT put your fingers in the mobile parts. The product does not contains any dangerous substance for the health. We disclaim all liability for any personal injury resulting from inappropriate installation or use of this seat.

Fiche produit mise à jour le : 10/01/2024